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# Fascia Bell Manual

Fascia Bell is a tool of Self-therapy which integrates cupping, TuiNa, Physiotherapy, fascial release and sports therapy skills, enhancing blood circulation of muscle and fascia to relieve pain.




Category	Petite	Classic	Max
<b>Suction:</b>	Small	Moderate	Max
<b>Application:</b>	Smaller areas: wrist, sides of neck and face	Larger areas: neck, shoulder, back, arms, legs...etc.	

## FDA Approved

- Recognized as the highest food safety standard by the World Health Organization
- Made by Food-contact silicone, which is harmless, massage oil friendly and can be disinfected with alcohol.
- Withstand temperatures up to 200°C

## Maintenance

- Clean with a wet wipe, alcohol or diluted household cleaning agent
- Prolonged soaking in detergent or using corrosive detergent is not recommended

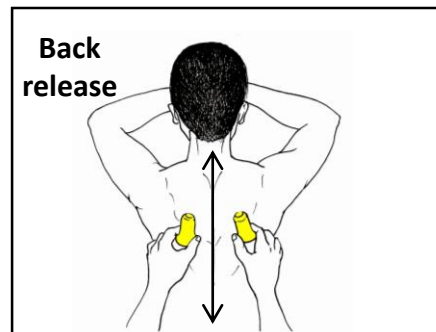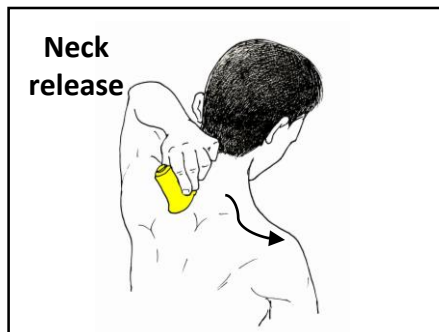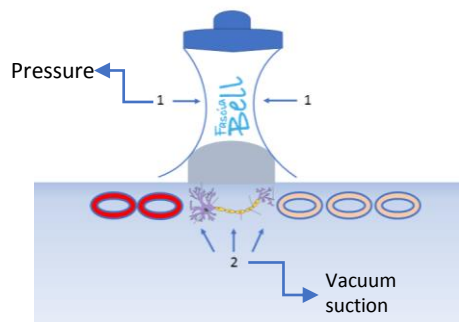Demonstration video:

Developed by:




Research Partner:

## Contraindication of using fascia bells

- Open wound, skin disease or allergy
- Acute pain or conditions of inflammation
- Recent fracture
- Illness with bleeding tendency
- Varicose veins
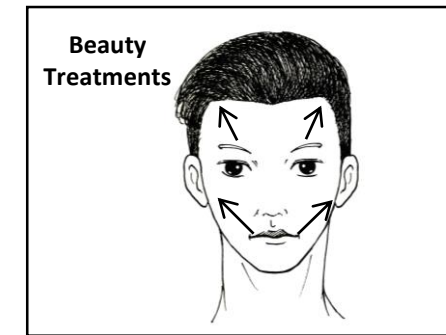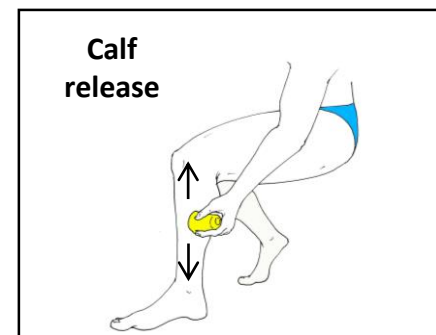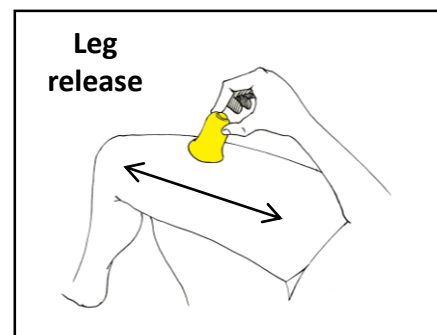
## Disclaimer

(i) Depend on individual condition, it is recommended to consult a health care professional when using this product. (ii) Ensure the Fascia Bell is intact before use. (iii) Keep yourself warm and avoid shower within 1-2 hours after using the Fascia Bell. (iv) In case of having blisters after using the Fascia Bell. Stay calm, disinfect the wound with alcohol and cover with gauze and immediately inform the professional for further treatment.